


CAP THEM ALL!


Now available! LVL fully-automated capper / decapper.
 The best choice for high throughput sample storage workflows.

- Available in SBS 96, 48 and 24 formats for all models of LVL SAFE® 2D barcoded tubes
- Automated height detection
- Capping/Decapping time respectively approx. 45 seconds

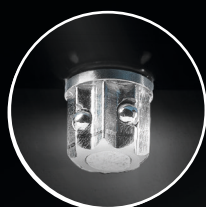

Latest cap driver technology: Cap drivers with spring-loaded metal balls optimally adapt to the caps and permit 100% leak-tight and material-saving sealing of the tubes without tilting.



 Technical specifications located on back

Find out more:
LVL cappers online



SAFE 
 Store Aliquots For Ever


Order now.

lvl-technologies.com/en/devices/capper--decapper

Innovative and low-maintenance rack-to-head technology:
It is not the complex driver head moves to the rack,
but the rack is lifted to the head.

- Can be operated as an integrated and standalone device
- Opens and closes a full rack of tubes at once
- Rack barcode recognition
- Touchscreen panel provides easy navigation through the system
- Caps can be hold or placed on a cap rack
- Ejects caps safely
- Very small footprint 408 x 234 x 388 mm

Optional accessory
LVL Cap Rack

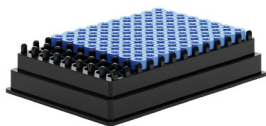

Conveniently store tube caps
on our cap racks.

SBS 96: RCK-CR96-P-L
SBS 48: RCK-CR48-P-L
SBS 24: RCK-CR24-P-L



96XT/96IT
fully-automated capper



48XT
fully-automated capper



24XT
fully-automated capper

Rack format	SAFE® SBS 96	SAFE® SBS 48	SAFE® SBS 24
Article number	External thread: CDC-96CH-XT-DD Internal thread: CDC-96CH-IT-DD	CDC-48CH-XT-DD	CDC-24CH-XT-DD